
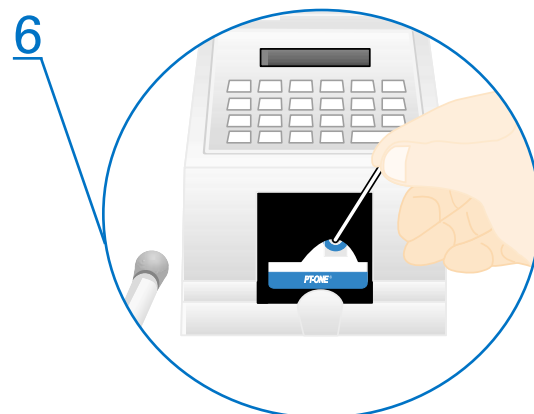


# Step-By-Step Procedures

1. Pass Card Through Reader
2. Insert Card
3. Input Your ID (Press ENTER)
4. Verify/Select Sample Type (Press ENTER)
5. Input Patient ID (Press ENTER)
6. Add Sample

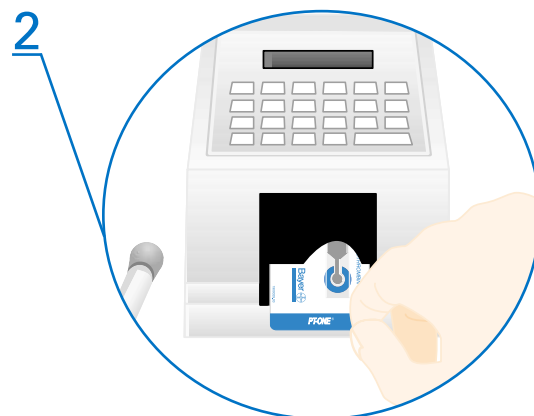

## Liquid Quality Control (LQC)

- Each test type has two levels
- Break ampule and mix vigorously 30 times
- Prepare and insert test card for QC
- Discard 5 drops into foil pouch
- Hold ampule vertically over testing area of card
- Apply one large drop to sample well
- Log QC results on appropriate log sheets
- QC results stored by test type

## Electronic Quality Control (EQC)

- 3-level test
- Each level has time and signal strength

| LEVEL | TIME (sec) | RANGE       | SIGNAL      |
|-------|------------|-------------|-------------|
| 1     | 30         | 29.1-30.9   | see package |
| 2     | 90         | 87.3-92.7   | insert for  |
| 3     | 270        | 261.9-278.1 | range       |

- To perform EQC, simply swipe card through reader and follow analyzer prompts
- No operator input needed for EQC ranges
- Results stored as 3 records



**NOTE** The manufacturer recommends running Liquid Quality Control once per week and Electronic Quality Control on every shift that performs coagulation testing.

**Part #**

118608  
118615  
119955  
120655  
120657  
120656  
119212  
123072  
128242

**Components**

Rapidpoint Coag Analyzer  
Rapidpoint Coag Power Supply/included with new unit  
EPSON® Printer Interface Cable  
EPSON Printer  
EPSON Printer Power Supply  
EPSON Printer Paper  
Rapidpoint Coag Instrument to PC Cable  
Rapidpoint ACCENT Analyzer  
Rapidpoint ACCENT Printer Paper

**Cards**

118564  
118570  
118580  
118590  
118599  
124366  
124367

PT-1.6ISI @ 100/box  
PT-1.0ISI @ 100/box  
PTNC-1.0ISI @ 100/box  
aPTT @ 100/box  
HMT @ 100/box  
PRT @ 25/box  
HTT @ 25/box

**Quality Control – Liquid & Electronic**

118576  
119151  
118586  
119152  
118595  
119153  
124430  
124431  
121907

RapidQC™ PT Level 1 @ 25/box  
RapidQC PT Level 2 @ 25/box  
RapidQC PTNC Level 1 @ 25/box  
RapidQC PTNC Level 2 @ 25/box  
RapidQC aPTT Level 1 @ 25/box  
RapidQC aPTT Level 2 @ 25/box  
RapidQC HMP Level 1 @ 25/box  
RapidQC HMP Level 2 @ 25/box  
EQC Test Card @ 1/box

**Miscellaneous & Service**

122016  
121039  
121038  
123236  
123237  
123234  
123235  
123232  
123233

TAS Upgrade for EQC  
Grey Reusable Pipettor  
Pipette Tips  
Demand Swap  
Demand Loan  
Annual Swap Contract  
Annual Loan Contract  
Extended Swap Contract  
Extended Loan Contract

**Contacts for Additional Assistance**

Technical Care Center: 1-877-BAYER11 (1-877-229-3711)  
Customer Service: 1-800-255-3232

Bayer Corporation  
Diagnostics Division  
511 Benedict Avenue  
Tarrytown, NY 10591

Rapidpoint is a registered trademark and RapidQC is a trademark of Bayer Corporation. EPSON is a registered trademark of Seiko Epson Corporation.  
©2002 Bayer Corporation. All rights reserved. Printed in USA. 0602.  
**0401242**


**Diagnostics**

**Making a positive  
difference to human health**

[www.bayerdiag.com](http://www.bayerdiag.com)